
Instruction manual & warranty for AIR BUBBLE Electric Cupping Massager

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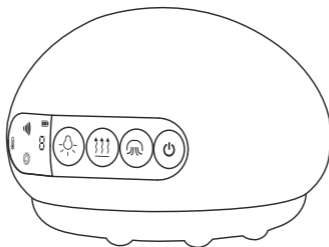www.habys.com

I. Specifications.

- **Product name:** AirBubble Electric Cupping Massager
- **Colour:** White
- **Voltage:** 3.7 V
- **Power:** 5 W
- **Heating levels:**
 - Low: 38 °C ± 3 °C (100.4 °F ± 5.4 °F)
 - Medium: 41 °C ± 3 °C (105.8 °F ± 5.4 °F)
 - High: 44 °C ± 3 °C (111.2 °F ± 5.4 °F)
- **Display:** LED
- **Intensity:** 5 levels
- **Noise level:** < 65 dB
- **Battery capacity:** 2200 mAh

II. Package Contents.

- 1 × AirBubble – Electric Cupping Massager
- 1 × USB-C cable
- 10 × Filter inserts
- 5 × Filter chambers
- 1 × User Manual and Warranty Terms


III. Features and Applications.

AirBubble – Electric Cupping Massager is a modern vacuum massage tool that combines traditional therapy techniques with innovative functions supporting relaxation, recovery, and body care. Thanks to its compact design, intuitive operation, and quiet performance, it is ideal for use in physiotherapy and beauty salons as well as at home.

Key features of the device:

- Vacuum massage with adjustable suction – supports lymphatic drainage, improves microcirculation, and facilitates toxin removal Heat therapy – three heating levels (38°C, 41°C, 44°C ±3°C) help relax muscles and promote recovery
- Red light therapy (infrared) – stimulates skin cells, enhances elasticity and firmness
- Magnetotherapy (8 magnetic points) – provides anti-inflammatory effects, improves circulation, and relieves tension
- Smart LED display – allows control of suction level, temperature, and functions
- Quiet operation (<65 dB) – ensures comfortable use in any environment
- Ergonomic design – lightweight, handy, and easy to use

IV. Device Description.

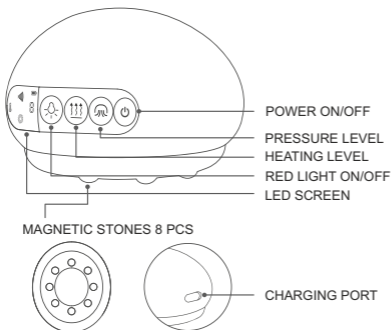V. Intended Use.

- Relaxation and recovery massage
- Reduction of muscle tension and fatigue
- Improved blood and lymph circulation
- Support for anti-cellulite care and body shaping
- Complement to cosmetic and physiotherapy treatments

VI. Safety Conditions for Use.

- Report any damaged cable or plug to the service center — do not attempt repairs yourself.
- Do not use the device while it is charging.
- Always disconnect the device from the power source before cleaning.
- Use only accessories recommended by the manufacturer.
- Do not pour liquids on the device — it is not waterproof.
- Avoid using the device in high-temperature or high-humidity environments (e.g., sauna, bathroom).
- Individuals with varicose veins, hypoglycemia, or hypertension should consult a physician before use.
- Do not use the product during menstruation or pregnancy.
- Wait 30–45 minutes after eating before using the device.
- After massage, wait at least 3 hours before bathing.
- Do not use the device daily — breaks between treatments are recommended.
- This product is not a medicine or medical device.

VII. Recommendations Before Using the Device.

- If you experience pain, discomfort, or skin irritation, stop the treatment and consult a specialist.
- Do not use the device on the chest, face, or neck.

- Avoid massaging the same area repeatedly — this may cause irritation or bruising.
- Perform the massage in one direction, following the flow of muscles or the lymphatic system.
- Remove all jewelry and other metal items before the session.
- Before use, apply a suitable oil to lightly moisturize the skin and ensure smooth gliding during the massage. If you have sensitive skin, perform an allergy test beforehand, e.g. on the wrist.

VIII. How to Use the Electric Cupping Device.

- **Turning on the device:** Press and hold the power button  until the red lights turn on — this indicates the device is ready for use.
- **Turning off the device:** Press and hold the power button  to switch off the device.
- **Adjusting pressure level:** A short press of the pressure button  changes the pressure level (1, 2, 3, 4, 5, 0 — in sequence). The pressure icon will flash, and a sound signal will confirm the change. A long press of the button  disables the pressure function.
- **Adjusting heating temperature:** A short press of the heating button  changes the heating level (1, 2, 3, 0 — in sequence). The heating icon will flash, and the change will be confirmed by a sound signal. A long press of the button  disables the heating function.
- **Red light function:** A short press of the red light button  turns the red light on or off.
- **Operating modes:**
 - Default mode: Red lights are illuminated.
 - Pressure mode: Levels 1–5.
 - Heating mode: Levels 1–3.
- **Heating temperature range:**
 - Low: 38°C ±3°C
 - Medium: 41°C ±3°C
 - High: 44°C ±3°C

IX. Electric Cupping Massage – Preparation and Techniques.

Before the massage, prepare the skin to ensure comfort and effectiveness. Apply an appropriate amount of massage oil to clean, dry skin. The oil facilitates smooth movement of the cup, protects against irritation, and enhances the vacuum effect. Always perform the massage in one direction — preferably upward, following the natural flow of muscles or the lymphatic system. Start with the lowest pressure level to allow the skin to gradually adapt and to assess its response.

Massage techniques for specific body areas:

- **Back:** Move upward, following the muscle lines.
- **Soles of the feet:** Massage from heel to toes with smooth, consistent strokes.
- **Abdomen:** Use circular motions in a clockwise direction or move upward.
- **Arms:** Massage from hands toward shoulders.
- **Legs, thighs, buttocks:** Begin at the feet or calves and move upward through the thighs to the buttocks and hips, following the muscle and lymphatic pathways.

X. Recommendations After Electric Cupping Massage.

- Avoid bathing, showering, or using a sauna for at least 3 hours — the skin needs time to regenerate.
- Clean the skin to remove any remaining oil and allow it to breathe freely.
- Drink a glass of water to support detoxification and recovery.
- Monitor your skin — if excessive redness, swelling, or pain occurs, discontinue use and consult a specialist.
- Do not perform the massage daily — the recommended frequency is 3–4 times per week.
- Keep the device in a dry place, clean the filter regularly, and never use the cup while charging.

XI. Charging the Device.

- The product is equipped with a lithium battery. When the display shows a low battery level, connect the device to a power source.
- Use a DC5V/1A adapter (not included) or a computer USB port.
- Connect the USB cable as follows:
 - Plug the USB end into the adapter or computer
 - Plug the USB-C end into the device
- The LED screen will display the charging status.
- The device must not be used while charging.
- When fully charged, the green indicator light stays on steadily.

XII. Cleaning and Maintenance.

- Before cleaning, make sure the device is turned off and unplugged from the power source.
- Wipe the casing with a soft, slightly damp cloth. For more stubborn dirt, you may use a disinfectant (max. 75%).
- Avoid any liquid contact with the internal components — this may cause a short circuit.
- The filter insert should be thoroughly rinsed or replaced with a new one, in accordance with the manufacturer's recommendations.

Filter Insert replacement:

1. Unscrew the filter cover to access the chamber where the filter is located.
2. Remove the used filter insert.
3. Insert the new filter element into the chamber.
4. Close the cover and ensure the filter is securely installed.

XIII. Storage.

Store the product in a location inaccessible to children, away from moisture, high temperatures, direct sunlight, and water splashes. If the device will not be used for an extended period, make sure to turn off the power.

XIV. Troubleshooting.

Issue	Possible Cause	Solution
Device does not turn on.	Battery is discharged.	Charge the device.
	Device is currently charging.	Use the device only after full charge.
Battery does not charge.	Plug is not properly connected.	Reconnect the plug.

XV. Battery – Usage and Disposal Guidelines.

The battery in this device is a dedicated rechargeable unit and is not suitable for use in other products.

Safety Guidelines:

- Keep away from fire and sources of high heat (e.g., oven, stove, direct sunlight).
- Do not disassemble or modify the battery.
- Avoid contact between the electrodes and metal objects.
- Do not store the device together with jewelry or other conductive materials.
- Do not charge or use the battery in overheated conditions.

XVI. Disposal.

- Devices marked with the crossed-out wheeled bin symbol must not be disposed of with household waste.
- They should be taken to a designated collection point for waste electrical and electronic equipment, in accordance with environmental protection regulations.

XVII. Legal Information.

- Compliance with EU regulations.
- The device meets the requirements of the following European Union directives:
 - Low Voltage Directive 2014/35/EU – EMC Directive 2014/30/EU
 - RoHS Directive 2011/65/EU
- The product is marked with the CE symbol.

XVIII. Warranty Terms:

1. For the purposes of this document:
 - a) "Manufacturer" refers to HABYS Sp. z o.o. (a limited liability company under Polish law), headquartered in Jasło, ul. Produkcyjna 16, 38-200 Jasło, registered in the National Court Register maintained by the District Court in Rzeszów, 12th Commercial Division, under KRS number: 0000513317, with a fully paid share capital of PLN 10,750,000 (ten million seven hundred fifty thousand zlotys), NIP: 6852208438, REGON: 180186290.
 - b) "Purchaser" refers to the end user of the equipment manufactured by the Manufacturer, including both consumers and business entities.
 - c) "Distributor" refers to an entity authorized by the Manufacturer to sell the Equipment within a designated territory.
 - d) "Authorized Service Center" refers exclusively to a service center operated by the Manufacturer or designated by the Manufacturer and operated by its Distributor.
 - e) "Equipment" refers to the movable item manufactured by the Manufacturer, subject to the warranty described herein.
2. **The Manufacturer guarantees the quality and proper functioning of the Equipment covered by this warranty card for a period of 2 years.**
3. Warranty repairs do not include activities specified in the user manual that the Purchaser is responsible for (e.g., assembly, maintenance).
4. Warranty liability covers only defects resulting from causes inherent in the Equipment. It does not cover damage due to improper use, storage, maintenance, or transport, nor discoloration due to sunlight exposure, chemical or abrasive agents, environmental contamination, accidental punctures or abrasions, stains from greasy or coloring substances (e.g., oils, creams), or permanent staining from dyes such as ink, permanent markers, or unstable textile dyes (e.g., blue jeans).
5. The warranty does not cover mechanical damage to the Equipment.
6. The Purchaser loses warranty rights if:

- a) the Equipment was repaired outside the Authorized Service Center during the warranty period,
 - b) any structural modifications were made without prior written consent from the Manufacturer,
 - c) the Equipment was used contrary to its intended purpose,
 - d) the Equipment was used after a structural defect was identified.
7. Regardless of the loss of warranty rights, the Manufacturer is not liable for any damage to property or harm to third parties resulting from use of the Equipment.
8. The warranty does not exclude or limit rights under statutory warranty (rękojmia). Exercising warranty rights suspends the statutory warranty period from the date the Manufacturer is notified of the defect until the Manufacturer refuses to fulfill warranty obligations or the deadline passes without resolution.
9. Warranty coverage is limited to the territory where the Equipment was sold.
10. The Manufacturer will fulfill warranty obligations (repair of defects caused by internal faults) free of charge within 14 days of receiving the Equipment. If repair cannot be completed within this period, the Manufacturer will notify the Purchaser and provide reasons and an estimated completion date.
11. To exercise warranty rights, the Purchaser must meet all of the following conditions:
- a) Provide a VAT invoice, receipt, or other proof of purchase, along with photographic documentation of the defect and a brief description of the claim.
 - b) Report the defect within 14 days of its discovery.
 - c) Deliver the Equipment in its original packaging at the Manufacturer's expense to the address listed in the warranty card or to the Authorized Service Center.
 - d) Indicate the return delivery address within the country of purchase.
12. The warranty period is extended by the duration of the repair, counted from the date the Equipment is delivered to the Manufacturer or Authorized Service Center until it is returned to the Purchaser.

13. The Manufacturer or Authorized Service Center determines the scope and method of repair.
14. The warranty is limited to repair. The Manufacturer reserves the right to replace the Equipment if the cost of repair exceeds the purchase price. The Purchaser is not entitled to demand replacement with a new item. The warranty does not cover losses or damages resulting from inability to use the Equipment.
15. The Manufacturer is not responsible for damage incurred during transport until the Equipment is received from the carrier. 16. If shipment was arranged by the Manufacturer, the Purchaser must inspect the Equipment upon delivery. If damage is found, the Purchaser must:
 - a) Notify the Manufacturer within 3 days of delivery to enable the Manufacturer to file a claim with the carrier.
 - b) Provide a damage report signed by the Purchaser and the carrier.
17. The warranty terms outlined in this document constitute the sole and exclusive warranty conditions for products of HABYS Sp. z o.o. (a limited liability company under Polish law).

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